



Example Protocol for the Culture of the Breast Carcinoma MDA-MB-231 Cell Line on Alvetex[®] Scaffold in Well Insert and Well Plate Formats

Introduction

MDA-MB-231 is a metastatic human breast cancer cell line originally isolated in the early 1970s[1]. MDA-MB-231 cells possess epithelial-like morphology and exhibit invasive properties when cultured *in vitro* and when transplanted into mice. MDA-MB-231 is widely used as a model for estrogen receptor negative breast cancer with applications in the study of tumorigenicity, metastasis, and cell invasion.

Alvetex Scaffold is currently available in four different cell culture formats: 24-well plate (AVP006), 12-well plate (AVP002), 6-well insert (AVP004), and 12-well insert (AVP005). 24-well and 12-well plates are suitable for shorter term cultures and for applications where limited cell penetration into the scaffold is required. Well insert formats generally support longer term cultures and deeper cell penetration into the scaffold. They also provide for conveniently tailored media set ups (see Quick Start Protocol). 6-well inserts can be placed in conventional 6-well plates, while 12-well inserts can be placed in either 6-well plates or 12-well plates, depending on media requirements. Alternatively, both insert types can be housed in the dedicated Well Insert Holder in Deep Well Petri Dish (AVP015) to allow for increased media volumes and prolonged cell culture. The availability of two different Alvetex Scaffold sizes enables choice on the basis of desired culture size and cell expenditure.

Methods

Preparation for 3D Cell Culture on Alvetex Scaffold

1. MDA-MB-231 cells (ATCC, HTB-26) were routinely maintained in T-75 flasks.

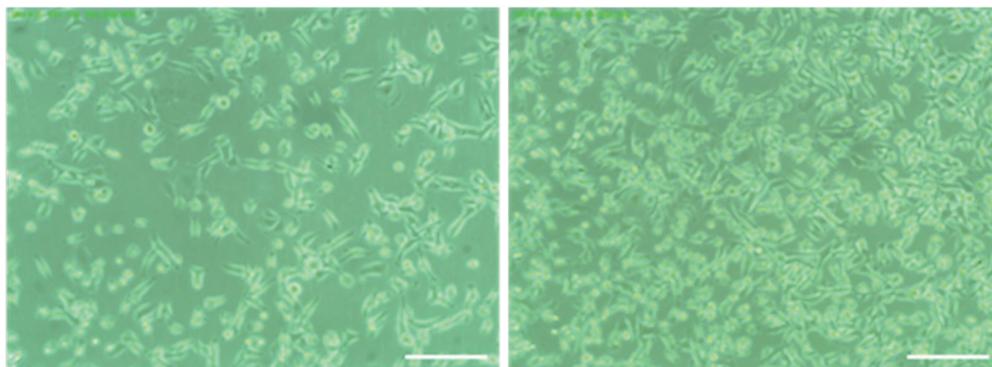


FIGURE 1. Phase contrast micrograph of MDA-MB-231 cells grown in conventional 2D culture plates. Image shows cells at low (left) and high (right) confluency. Scale bars: 100 μ m.

2. Complete growth media consisted of: Leibovitz's L-15 medium (ATCC, 30-2008) supplemented with 10 % v/v FBS and 100 U/ml Penicillin/Streptomycin.

Note: The L-15 medium already contains 2 mM L-glutamine and is formulated for use in air of atmospheric composition.



Example Protocol for the Culture of the Breast Carcinoma MDA-MB-231 Cell Line on Alvetex[®] Scaffold in Well Insert and Well Plate Formats

3. Cells were harvested by trypsinisation and centrifuged for 5 minutes (1000 rpm). The supernatant was discarded and the cell pellet was re-suspended in an appropriate volume of media for cell counting by Trypan Blue.
4. Cells were re-suspended at a concentration of 3.33×10^6 cells/ml for seeding.

24-well Plate Format, AVP006

1. Alvetex Scaffold 24-well plates were prepared for seeding with a 70 % ethanol wash (2 ml per well) and subsequent media washes (twice with 2 ml of media each).
2. 75 μ l of the cell suspension was added to the centre of each Alvetex Scaffold disc, which was equivalent to 0.25×10^6 cells per well.
3. The plate was incubated for at least 15 minutes at 37 °C without CO₂ to allow the cells to settle into the scaffold.
4. 2 ml of media was added to each well taking care not to dislodge cells from Alvetex Scaffold.
5. Plates were re-incubated and maintained by complete media exchange every other day.

Note: This method can be applied to the use of Alvetex Scaffold in 12-well plate format, AVP002. Adjust cell seeding and media volumes according to the guidelines provided in the Quick Start Protocol.

12-well Insert Format, AVP005

1. Alvetex Scaffold 12-well inserts in 12-well plate format were prepared for seeding by dipping in 70 % ethanol followed by media washes (twice with 4 ml per well).
2. 75 μ l of the cell suspension was added to the centre of each Alvetex Scaffold disc, which was equivalent to 0.25×10^6 cells per well.
3. The plate was incubated for at least 15 minutes at 37 °C without CO₂ to allow the cells to settle into the scaffold.
4. Media was carefully added to the outer wall of each well to a total volume of 4 ml, taking care not to dislodge cells from Alvetex Scaffold.
5. Plates were re-incubated and maintained by complete media exchange every other day.

Note: This method can be applied to the use of Alvetex Scaffold in 6-well insert format, AVP004. Adjust cell seeding and media volumes according to the guidelines provided in the Quick Start Protocol.



Example Protocol for the Culture of the Breast Carcinoma MDA-MB-231 Cell Line on Alvetex[®] Scaffold in Well Insert and Well Plate Formats

Example Data: 24-well plate format, AVP006

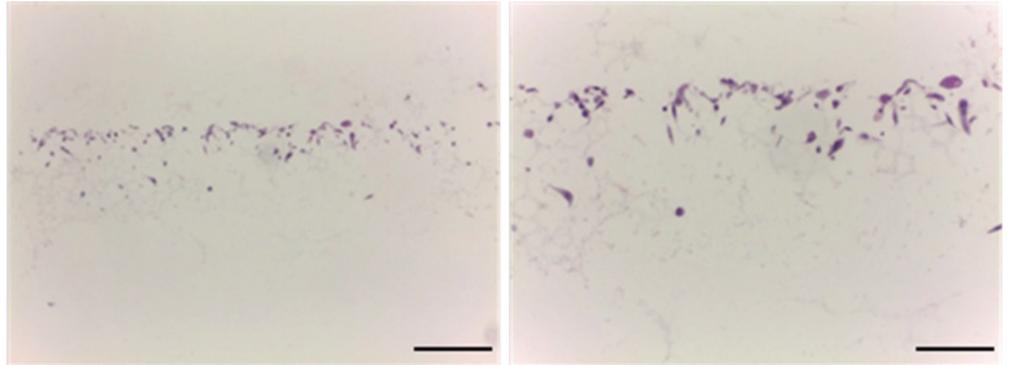


FIGURE 2. Brightfield micrograph showing the structure of MDA-MB-231 cells cultured for 3 days on 15 mm diameter Alvetex Scaffold discs presented in 24-well plate format. Cells were fixed, embedded in paraffin wax, sectioned (10 µm) and counterstained with haematoxylin and eosin. Scale bars: Left: 200 µm, Right: 100 µm.

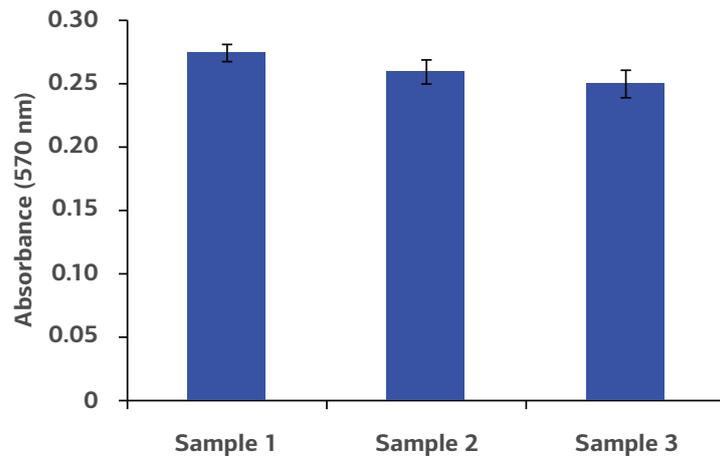


FIGURE 3. Biochemical analysis of cell viability using a standard MTT assay. Data from 2 sample replicates of MDA-MB-231 cells are shown, each sampled in triplicate (n=3, mean ± SD). Cells were cultured for 3 days on 15 mm Alvetex Scaffold discs presented in 24-well plate format.



Example Protocol for the Culture of the Breast Carcinoma MDA-MB-231 Cell Line on Alvetex[®] Scaffold in Well Insert and Well Plate Formats

Example Data: 12-well insert format, AVP005

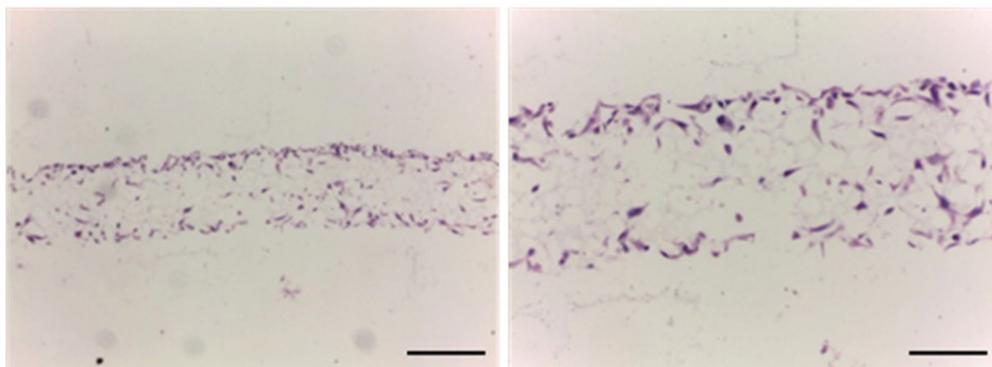


FIGURE 4. Brightfield micrograph showing the structure of MDA-MB-231 cells cultured for 3 days on 15 mm diameter Alvetex Scaffold discs presented in 12-well insert in 12-well plate format. Cells were fixed, embedded in paraffin wax, sectioned (10 μ m) and counterstained with haematoxylin and eosin. Scale bars: Left: 200 μ m, Right: 100 μ m.

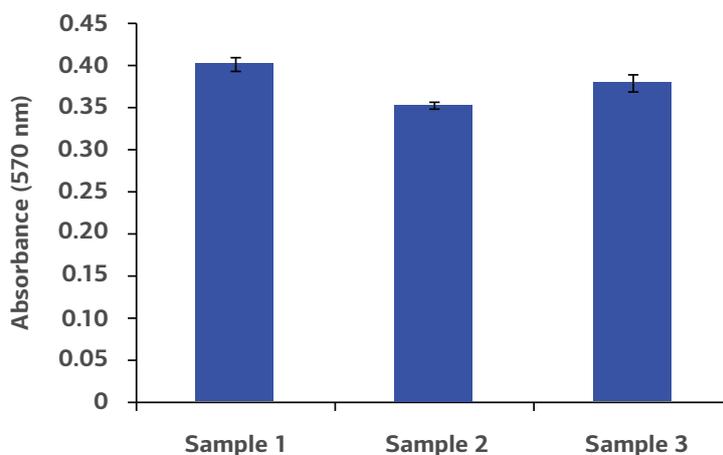


FIGURE 5. Biochemical analysis of cell viability using a standard MTT assay. Data from 2 sample replicates of MDA-MB-231 cells are shown, each sampled in triplicate (n=3, mean \pm SD). Cells were cultured for 3 days on 15 mm Alvetex Scaffold discs presented in 12-well inserts in 12-well plate format.

References

1. Cailleau, R.M. *et al.*, 1974. Breast tumor cell lines from pleural effusions. *J Natl Cancer Inst* 53 (3): 661-74.



Example Protocol for the Culture of the Breast Carcinoma MDA-MB-231 Cell Line on Alvetex[®] Scaffold in Well Insert and Well Plate Formats

Technical Support

To obtain the latest technical support information, please visit www.reinnervate.com. You can have 24-hour access to the following information:

- Access to telephone and fax information for order support
- Search through the Reinnervate protocol video gallery
- Search product spec sheets, protocols, applications notes and product support docs
- Obtain information on customer training and workshops
- Submit a question directly to Technical Support at tech.support@reinnervate.com

Ordering Information

For the latest products and services available through Reinnervate, please visit www.reinnervate.com.

Reinnervate

NETPark Incubator, Thomas Wright Way,
Sedgefield, Co. Durham, TS21 3FD
UNITED KINGDOM

T: +44 (0) 1740 625 266

F: +44 (0) 1740 625 251

E: info@reinnervate.com

Trademarks

Unless otherwise noted, all trademarks are the property of Reinnervate Ltd.
©2015 Reinnervate Ltd. All rights reserved.